
1C56 Taleo Adjust

DE	Gebruiksaanweisung	7
EN	Instructions for use	16
FR	Instructions d'utilisation	24
IT	Istruzioni per l'uso	34
ES	Instrucciones de uso	43
PT	Instruções de utilização	52
NL	Gebruiksaanwijzing	61
SV	Bruksanvisning	70
DA	Brugsanvisning	78
NO	Bruksanvisning	87
FI	Käyttöohje	95
PL	Instrukcja użytkowania	104
HU	Használati útmutató	113
CS	Návod k použití	121
RO	Instrucțiuni de utilizare	130
HR	Upute za uporabu	139
SL	Navodila za uporabo	147
SK	Návod na použitie	156
BG	Инструкция за употреба	164
TR	Kullanım kılavuzu	174
EL	Οδηγίες χρήσης	183
RU	Руководство по применению	193
UK	Інструкція щодо використання	202
JA	取扱説明書	211
ZH	使用说明书	219
KO	사용 설명서	226

Bench alignment, transtibial

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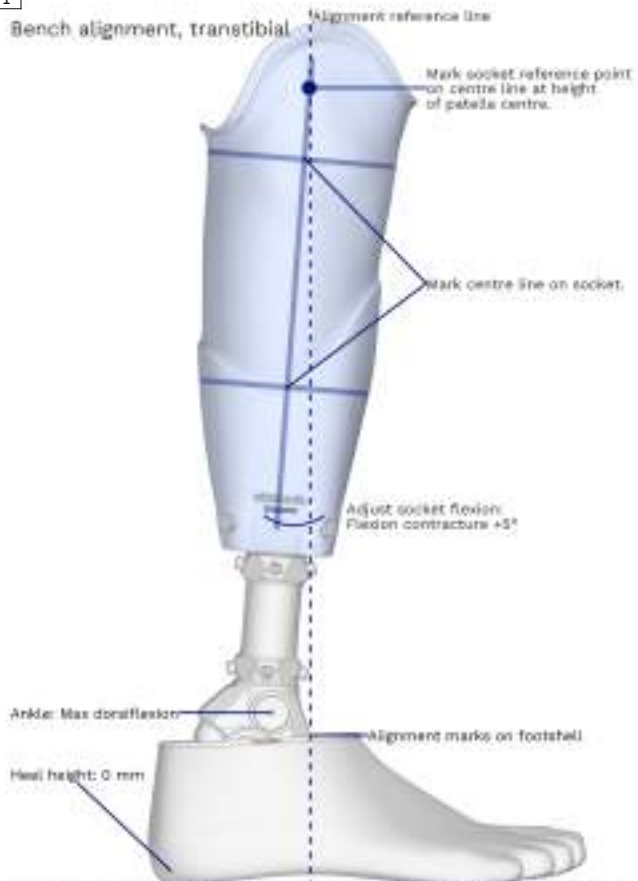INFORMATION

Date of last update: 2025-03-07

- ▶ Please read this document carefully before using the product and observe the safety notices.
- ▶ Instruct the user in the safe use of the product.
- ▶ Please contact the manufacturer if you have questions about the product or in case of problems.
- ▶ Report each serious incident related to the product to the manufacturer and to the relevant authority in your country. This is particularly important when there is a decline in the health state.
- ▶ Please keep this document for your records.

The heel height of the Taleo Adjust prosthetic foot 1C56 is infinitely adjustable (**0 mm** to **70 mm**) to enable shoes with different heel heights to be worn.

Carbon and polymer spring elements permit perceptible plantar flexion at heel strike, a natural rollover and high energy return. A removable heel wedge absorbs the heel strike.

2 Safety-related information

CAUTION! Risk of injury and risk of product damage

- ▶ To prevent mechanical damage, treat the product with care.
- ▶ Observe the information on product combinations in the instructions for use for the products.
- ▶ Observe the maximum lifetime of the product (see page 18).
- ▶ Check the product for damage and readiness for use prior to each use.
- ▶ Observe the operating conditions of the product (see page 17) and do not expose the product to prohibited environmental conditions (see page 17).
- ▶ Do not use the product if it is damaged or in a questionable condition. Clean the product if it gets dirty. If necessary, have it inspected, repaired or replaced by the manufacturer or a specialist workshop.
- ▶ The product may only be used by one person; reuse on another person is prohibited.

Signs of changes in or loss of functionality during use

Reduced spring effect (e.g. decreased forefoot resistance or changed rollover) or delamination of the spring are indications of loss of functionality. Play in the hydraulic ankle joint, oil leakage and unusual noises may indicate a loss of functionality.

3 Intended use

3.1 Indications for use

The product is intended exclusively for lower limb exoprosthetic fittings.

3.2 Operating conditions

The MOBIS classification describes the mobility grade and body weight, and makes it easy to identify compatible components.

The product is recommended for mobility grade 2 (restricted outdoor walker) and mobility grade 3 (unrestricted outdoor walker).

- **Product not approved for:** Use in hip disarticulation prostheses

Stiffness

The table below shows the appropriate spring stiffnesses for the user's body weight.

Within the weight ranges, the user can choose between two spring stiffnesses according to preferences.

TIP: Users requiring extra stability can test the stiffness for a higher load level.

Body weight [kg]	Spring stiffness
up to 51	1 or 2
52 to 58	2 or 3
59 to 67	3 or 4
68 to 77	4 or 5
78 to 88	5 or 6
89 to 100	6 or 7
101 to 115	7

3.3 Environmental conditions

Allowable environmental conditions
Temperature range: -10 °C to +45 °C
Chemicals/liquids: Fresh water, soap solution, chlorine water, salt water ≤ 3.5 % salt content

Allowable environmental conditions
Moisture: Submersion: max. 1 h in 2 m depth, relative humidity: no restrictions
Solids: dust, occasional contact with sand
Clean the product after contact with humidity/chemicals/solids, in order to avoid increased wear and tear and damage (see page 23).
Prohibited environmental conditions
Solids: highly liquid-binding particles (e.g. talcum), dust in high concentrations (e.g. construction site), intensive contact with sand
Chemicals/liquids: perspiration, urine, acids, continuous use in liquid media
Storage and transport
Storage temperature: -20 °C to +60 °C, relative humidity: 20 % to 90 %, no mechanical vibrations or impacts

3.4 Product combinations

This prosthetic component is compatible with Ottobock's system of modular connectors. Functionality with components of other manufacturers that have compatible modular connectors has not been tested.

3.5 Lifetime

Prosthetic foot

Depending on the user's activity level, the maximum lifetime of the product is 3 years.

Footshell, protective sock

The product is a wear part, which means it is subject to normal wear and tear.

4 Scope of delivery

Quantity	Designation	Reference number
1	Instructions for use	–
1	Prosthetic foot	–
1	Protective sock	SL=Spectra-Sock-7
1	Heel wedge set	2F50=*

Spare parts/accessories	
Designation	Reference number
Footshell	2C8=*
Connection cover for 2C8=*	2C21=*
Footshell, narrow	2C8=*S

5 Preparing the product for use

CAUTION

Incorrect alignment, assembly or adjustment

Injury due to incorrectly installed or adjusted as well as damaged prosthetic components

- Observe the alignment, assembly and adjustment instructions.

CAUTION

Reaching into the ankle joint

Risk of injury due to pinching of the fingers

- Do not reach into the ankle joint area.

NOTICE

Altering the prosthetic foot or footshell

Premature wear resulting from damage to the product

- Do not alter the prosthetic foot or footshell.

INFORMATION

A plastic part protects the connecting section of the product from scratches during the alignment and fitting.

- Remove the plastic part before the patient leaves the fitting area.

5.1 Applying/removing the footshell

- Pull the protective sock over the prosthetic foot to prevent noises in the footshell.
- **CAUTION! Always use the prosthetic foot with the footshell.**
Apply or remove the footshell as described in the footshell instructions for use.

5.2 Bench Alignment

INFORMATION: Use the prosthetic foot only with an installed heel wedge.

INFORMATION: If the user wants more stability: Move the prosthetic foot approx. 1 cm further in the anterior direction. This lengthens the forefoot and increases the knee-securing effect.

- > **Required materials:** goniometer 662M4, heel height measuring device 743S12, 50:50 gauge 743A80, alignment apparatus (e. g. PROS.A. Assembly 743A200)
 - > **Illustration to prepare for bench alignment (see fig. 1)**
 - > **Transfemoral prosthesis or knee disarticulation prosthesis:** Observe the instructions for use for the prosthetic knee joint.
- 1) Position the prosthetic foot in the alignment apparatus.
 - **Heel height: 0 mm** (barefoot)
 - **Ankle joint:** Bring the prosthetic foot into maximum dorsiflexion (see page 22) so the complete adjustment range can be accessed.
 - 2) **Sagittal plane:** Align the prosthetic foot with the alignment reference line.
 - Alignment reference line through the markings of the footshell.
 - 3) **Frontal plane:** Align the prosthetic foot with the alignment reference line.
 - Alignment reference line through the markings of the footshell (exterior rotation **approx. 5°**, alignment reference line in the middle of the second toe).
 - 4) Install the prosthetic components between the prosthetic foot and prosthetic socket.
 - 5) **Sagittal plane:** Align the prosthetic socket with the alignment reference line.
 - Determine the centre of the prosthetic socket proximally and distally with the 50:50 gauge and mark the centre line.
 - Mark the socket reference point on the centre line: **For transtibial prostheses, at the height of the centre of the patella.**
 - Align the prosthetic socket so the alignment reference line runs through the socket reference point.
 - Adjust the socket flexion by rotating around the socket reference point: **individual residual limb flexion/flexion contracture + 5°**
 - 6) **Frontal plane:** Align the prosthetic socket with the alignment reference line.
 - Take the abduction or adduction of the residual limb into account.

5.3 Static Alignment

The static alignment of the Taleo Adjust consists of two stages. You have to perform both stages in order to set the minimum and maximum heel height.

Ottobock recommends checking the alignment of the prosthesis during both stages using the L.A.S.A.R. Posture.

1) Without footwear (see fig. 2):

- To perform the static alignment without footwear, bring the prosthetic foot into maximum dorsiflexion.
- Use the set screws to make any necessary adjustments.

2) With everyday footwear (see fig. 3):

- Carry out the static alignment with the everyday footwear.
- Only adjust the angle of the sagittal plane by changing the heel height (see page 22), not using the set screws. Using the set screws changes the maximum/minimum adjustable heel height.

The alignment recommendations (TF modular leg prostheses: **646F219***, TT modular leg prostheses: **646F336***) can be requested from Ottobock as needed.

5.4 Dynamic Trial Fitting

The optimal gait pattern is developed during dynamic fitting. The prosthesis is aligned optimally in the frontal plane and the sagittal plane for this purpose.

- ▶ Use everyday footwear for dynamic fitting.
- ▶ **TT fittings:** Make sure that physiological knee movement in the sagittal and frontal planes is achieved after the leg begins to bear weight following the heel strike. Avoid medial movement of the knee joint.
 - If the knee joint moves in the medial direction in the first half of the stance phase, move the prosthetic foot in the medial direction.
 - If the knee joint moves in the medial direction in the second half of the stance phase, reduce the exterior rotation of the prosthetic foot.
- ▶ Once the gait pattern has been optimised with everyday footwear, try additional heel heights. The heel height affects rollover as it changes the angle at which the prosthetic foot is loaded.
- ▶ After completing the dynamic fitting and the walking exercises: Remove the plastic part from the pyramid.

5.4.1 Replacing the heel wedge

The behaviour of the prosthetic foot at heel strike and during heel contact in the mid-stance phase can be adapted by replacing the heel wedge. Heel wedges of various hardness are included in the scope of delivery.

Heel wedge hardness: The colour of the heel wedge indicates the hardness (see fig. 4). Ottobock recommends starting with the preinstalled heel wedge.

- 1) Bend the heel spring of the prosthetic foot slightly and remove the installed heel wedge.
- 2) Align the other heel wedge so the Ottobock lettering is right side up and the tip points in the anterior direction.
- 3) Insert the heel wedge into the prosthetic foot (see fig. 7). To do so, first insert the tip into the front part of the retainer ❶, then push the heel wedge into the rear of the retainer ❷.

6 Use

6.1 Setting the heel height

INFORMATION

First-time use or after an extended break: there may be stiffness when operating the lock axle and/or adjusting the ankle angle.

- To loosen: slide the lock axle back and forth several times and move the ankle joint through its entire scope of movement.

The heel height of the prosthetic foot can be smoothly adjusted from **0 mm** to **70 mm**.

Adjustment method

Set the heel height while seated. Then hold on to a suitable object, stand up and test whether the set heel height is appropriate or needs to be corrected again.

Setting the heel height

- 1) Note how your prosthesis feels when you walk barefoot.
- 2) Sit down.
- 3) Note the angle of your prosthetic leg when you sit barefoot and the prosthetic foot is flat on the floor.
- 4) Put on the shoes you want to wear now.
- 5) **CAUTION! Do not put any weight on the prosthesis while setting the heel height. There is a risk of injury given the freely movable nature of the prosthetic foot.**
Unlocking the ankle joint: Push the lock axle in from the inside of the ankle as far as it will go (see fig. 8). The ankle joint can now be moved.
- 6) Push the tip of the prosthetic foot towards the sole of the foot as far as it will go (plantar flexion).

- 7) Place your prosthetic leg on the floor. The shoe should be flat on the floor and the prosthetic leg should have the same angle as it has when you are sitting barefoot.
- 8) **Locking the ankle joint:** Push in the lock axle protruding from the outside of the ankle as far as it will go (see fig. 9). The ankle joint is now locked and can bear weight once more.
- 9) Hold on to a suitable object, stand up and test whether the heel height setting is appropriate.

7 Cleaning

- > **Allowable cleaning agent:** pH neutral soap (e.g. 453H10 Derma Clean)
- 1) **NOTICE! To avoid product damage, only use the permissible cleaning agents.**
Clean the prosthetic foot and footshell with clear fresh water and pH-neutral soap.
- 2) **If present:** Use a toothpick to remove dirt from the contours that drain off the adapter connections or the sole of the footshell and rinse.
- 3) Rinse the soap away with clear fresh water. In doing so, rinse the footshell until all dirt has been removed.
- 4) **In case of heavy soiling and occurring sounds:** Pull the footshell and protective sock off the prosthetic foot. Clean all products individually with water.
- 5) Dry the product with a soft cloth.
- 6) Allow to air dry in order to remove residual moisture.

8 Maintenance

- ▶ A visual inspection and functional test of the prosthetic components should be performed after the first 30 days of use.
- ▶ During regular inspection: Inspect the prosthesis for signs of wear and check its function.
- ▶ Conduct annual safety inspections.

Information on the sacrificial anode

There is a sacrificial anode made of zinc in the pyramid that increases the corrosion protection of the product. The sacrificial anode is maintenance-free and can fully wear over the lifetime of the product.

9 Disposal

Do not dispose of the product in household waste. Improper disposal is harmful to health and the environment. Observe the instructions of the local authority regarding return, collection and disposal.

10 Legal information

All legal conditions are subject to the respective national laws of the country of use and may vary accordingly.

10.1 Liability

The manufacturer will only assume liability if the product is used in accordance with the descriptions and instructions provided in this document. The manufacturer will not assume liability for damage caused by disregarding the information in this document, particularly due to improper use or unauthorised modification of the product.

10.2 CE conformity

The product meets the requirements of Regulation (EU) 2017/745 on medical devices. The CE declaration of conformity can be downloaded from the manufacturer's website.

11 Technical data

Sizes [cm]	22	23	24	25	26	27	28
Heel height [mm]	0 to 70						
System height with 2C8*, 2C8*S [mm]	92		95		100		
Build height with 2C8*, 2C8*S [mm]	110		113		118		
Weight with 2C8* [g]	710	750	770	815	855	880	925
Weight with 2C8*S [g]	705	740	760	805	840	860	–
Max. body weight [kg]	88	100		115			
Mobility grade	2, 3						

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